MANUALE D'USO SEDIA DOCCIA AQUA ATLANTIC

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USER MANUAL SHOWER CHAIR AQUA ATLANTIC

INTRODUCTION

Thank you for choosing Demarta Virginio. The products have been designed and manufactured to meet all your needs for practical, correct and safe use. This manual contains small suggestions for the correct use of the device you have chosen and valuable advice for your safety. We recommend that you read this manual carefully before using the item you purchased. If in doubt, please contact the retailer, who will be able to help you and advise you correctly. The manufacturer reserves the right to make changes to the device and the following manual without notice in order to improve its characteristics.

GENERAL WARNINGS

Use the product in the manner indicated in this manual. If you still have doubts, contact DEMARTA VIRGINIO S.r.l.

Always evaluate the compatibility of the product with the characteristics of the subject together with your doctor or therapist.

Always check the integrity of the product before each use, even when new. Do not use if damaged or worn.

DEMARTA VIRGINIO S.r.l. declines all responsibility for damage to property or persons resulting from improper use, not agreed with a professional, not corresponding to what is indicated in the user manual.

Any accidents, anomalies, breakages must be reported to DEMARTA VIRGINIO S.r.l. in order to monitor the quality and critical issues of the product.

Keep the packaged product away from any heat source as the packaging is made of cardboard.

Always pay attention to the presence of children.

The product has been designed only for the intended use. Do not abuse it in any way.

EC DECLARATION OF CONFORMITY

DEMARTA VIRGINIO S.r.l. declares under its sole responsibility that the product manufactured and placed on the market by DEMARTA VIRGINIO S.r.l. complies with the applicable provisions of EU Directive 2017/745 on medical devices.

To this end, DEMARTA VIRGINIO S.r.l. guarantees and declares under its sole responsibility the following:

The devices meet the general safety and performance requirements as required by Annex I of Regulation 2017/745 as prescribed by Annex IV of the aforementioned regulation.

- The devices ARE NOT MEASURING INSTRUMENTS.
- The devices ARE NOT INTENDED FOR CLINICAL INVESTIGATIONS.
- The devices are marketed in NON-STERILE packaging.
- The devices are to be considered as belonging to class I in compliance with the provisions of Annex VIII of the aforementioned regulation.

INTENDED USE

Evacuation and mobility aid intended for people with walking difficulties in the bathroom.

SPECIFIC WARNINGS

Always be very aware of moving parts that could cause limb entrapment and personal injury.

Always consult your doctor or therapist for proper use of the device.

The service life of the device is determined by wear and tear of non-repairable and/or replaceable parts.


TECHNICAL FEATURES

Width (cm)	Seat Width (cm)	Seat Height (cm)	Maximum Load (Kg)	Net Weight (Kg)	Carton Dimensions (cm)
					
50	43	49 / 58	120	11.5	52 x 20 x 54

ASSEMBLY

1) Frame assembly.

Remove the white tube from the packaging and insert it into the hole provided at the two ends of the frame.

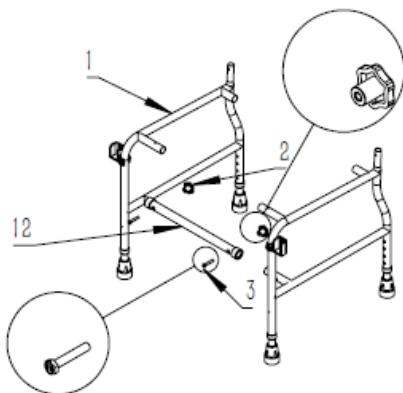After inserting the tube, secure the frame by inserting the 2 star locking knobs into the holes.

ATTENTION!


After attaching the 2 star knobs, make sure of the correct stability of the frame.

2) Backrest assembly.

Remove the backrest shown in the image below from the box and insert it into the holes provided.

Fix the two backrest tubes to the frame by inserting the two clips supplied into the appropriate holes.

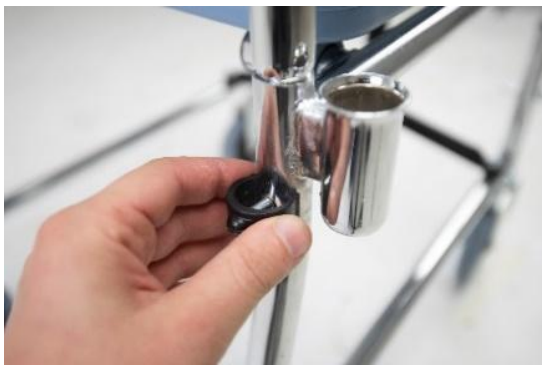3) Armrest assembly.

Remove the two side armrests from the packaging.

Insert the two side armrests into the appropriate slots on the outside of the seat.

ATTENTION!

Be sure to insert the two undrilled tubes near the backrest first, then attach the other two into the front of the shower chair.

To facilitate the transport of the chair and the transfer of the patient, all armrests are completely foldable.

ATTENTION!

Never use the armrests to raise or move the chair.

4) Seat assembly.

Remove the seat from the packaging.

Place it on top of the frame and apply light pressure on it until the 4 clips are correctly locked.

ATTENTION!

Make sure it is properly secured before proceeding with assembly.


5) Inserting the seat end cap.

Remove the seat cap from the packaging and insert it into the hole provided between the basin and the seat.

6) Tray assembly.

Remove the tray from the packaging and place it under the seat, Insert it by sliding it from the back side of the chair, taking care to insert it correctly into the appropriate guides.

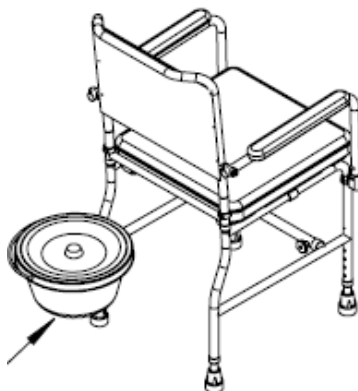7) Leg height adjustment.

The legs of the shower chair have several holes where you can adjust the height.

- Press the metal spring button so that it retracts fully into the tube.
- Select the desired height by adjusting the handles in an upward or downward motion.
- Release the metal spring button in front of the desired hole location.

ATTENTION!

Before use, always make sure that each spring-loaded button protrudes completely from the selected hole position and that the height is equal for both 4 legs.

BEFORE EACH USE

- Always check the state of wear and tear of the mechanical parts involved so as to ensure that the device can be used in total safety for people and property.
- Check that the device is assembled correctly.

WARNINGS FOR USE

- Do not exceed the maximum permitted capacity (120 kg) for any reason.
- Always assist the sick person during handling operations.
- Before transferring the patient to/from the chair, make sure that the chair is placed on a flat, stable surface with the wheels locked.
- To avoid possible tipping, do not lean over the back of the chair.
- Keep the chair away from direct heat sources.

USE

Getting on the shower chair

- Bring the chair as close as possible so that it is positioned next to you.
- Lean against the frame structure by forcing yourself on your arms, then slowly descend towards the seat.

Getting off the shower chair

- Get as close as possible to a new seat.
- Place your arms on the frame of the seat on which you are to move and slowly slide towards the latter

ATTENTION!

Assistance or supervision by experienced and qualified personnel is required to prevent slipping and falls.

MAINTENANCE

The shower chair requires a few elementary periodic maintenance, below we list the checks to be carried out periodically.

- Cleaning of parts (seat, backrest, armrests, wheels, etc.) as needed.
- Pin lubrications: monthly.

It is of fundamental importance that the shower chair is completely checked in all its parts when used for the first time, in particular that the screws and nuts are all tightened.

CLEANING AND DISINFECTION

Cleanliness

Use a cloth soaked in warm water and mild soap and dry with a clean, dry cloth.

Disinfection

If necessary, disinfect the product using a mild disinfectant detergent for normal household use. The toilet bowl should be regularly cleaned and disinfected using a normal bath disinfectant.

ATTENTION!

Do not use abrasives, acids, alcohol, chlorine-based detergents, disinfectants and acetone, as these substances cause abrasion of the shower chair components.

The manufacturer will not be responsible for damage caused by the use of materials that may deteriorate the surface of the product or by chemical corrosives used during cleaning.

GENERAL DISPOSAL CONDITIONS

When disposing of the device, never use normal municipal solid waste delivery systems. It is recommended to dispose of the device through the common municipal ecological islands for the planned recycling operations of the materials used.

WARRANTY

All Demarta products are warranted against defects in material or workmanship for a period of 2 (two) years from the date of sale of the product, subject to any exclusions and limitations specified below. This warranty is invalid in the event of misuse, abuse or modification of the product and for failure to follow the operating instructions. Demarta does not warrant the products against damage or defects in the following conditions: natural disasters, failure to follow the guidelines and instructions for use, or failure to perform maintenance as indicated in the manual.

Components that can wear and tear are not covered by this warranty if the damage is due to normal use of the product.